


# Quick Start Guide

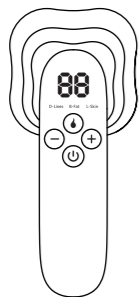

## Cellulite Massager F2

Iron off stubborn fat and stretch marks

- Thank you for using the cellulite massager F2. Please read Quick Start Guide before using, and follow its instruction.
- In order to use this product safely, please read the "Important Safety Tips" in the Quick Start Guide carefully.
- After reading, please keep manual in a safe place for reference at any time.
- Please keep the product out of the reach of children.
- This product belongs a non-medical device.

### Important Safety Tips

Tips	
Users Disabled	<ul style="list-style-type: none"> <li>• Patients who use pacemakers or have serious arrhythmias.</li> <li>• Patients with cognitive impairment and skin sensory paralysis.</li> <li>• Babies, pregnant women, parturients (60 days after surgery).</li> </ul>
Parts Forbid -en	<ul style="list-style-type: none"> <li>• Implanted sites such as metal and plastic.</li> <li>• Skin ulcers.</li> <li>• The heart, both carotid sinuses and the front of the throat.</li> </ul>
Action Prohibit -ed	<ul style="list-style-type: none"> <li>• Do not immerse this product in water.</li> <li>• It is forbidden to use it while operating the other machine or driving.</li> <li>• It is prohibited to disassemble, modify and repair this product by yourself.</li> </ul>
Instructi -ons Followed	<ul style="list-style-type: none"> <li>• This product should not be near strong magnets.</li> <li>• Do not use it on the same part of the body for a long time, it may cause skin inflammation.</li> <li>• The charging cable cannot be replaced by yourself. If it is damaged, please contact the customer for confirmation before replacing it.</li> </ul>

### Instructions For Use


Please follow the steps below as first-time usage:

1. Remove the protective film on the electrode sheet and handle on the device.
2. Wipe off the remaining glue on the electrode pad and the protective film on the handle with a wet paper towel firstly, then wipe off the water marks with a dry paper towel.
3. Power on and confirm whether each function key works normally.

### Tips:


Please do not dispose of the product randomly after the end of its service life. The built-in battery may contain environmentally harmful ingredient. Please dispose of the product at local disposal bins and make sure that the internal power of the battery is exhausted.

### Product Specification

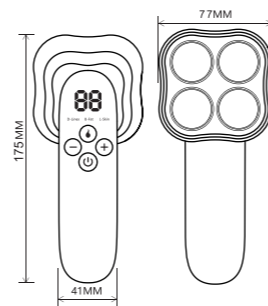



Battery Capacity	3.7V/1200mAh	RF Frequency	1M HZ
Charge	5V/2A	Red Wavelengt	630(±10)nm
BIO-Wave Frequency	30-300HZ	Material	ABS, PC, Stainless Steel
Temperature Control	37-42℃	Package Weight	265g
Product Weight	130g		

### Accessories Details

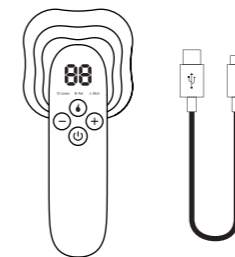

Device	USB charging cable
1	1

### Component Function

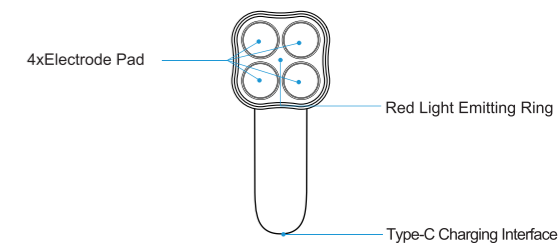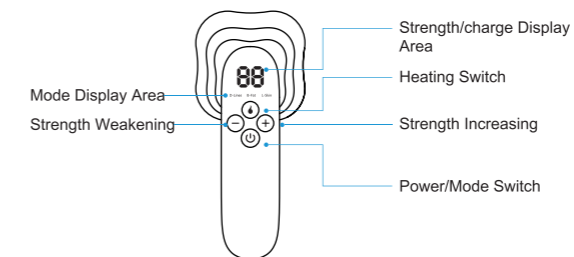

### Charge

- **PMU display:**
  - Charging: ring dynamic display.
  - Charging overvoltage: display OC, please replace the adapter compatible with for charging.
  - Fully charged: display FU.
  - Low voltage: display LO.
- Insert the Type-C connector that comes with the device into the charging port of the device, then plug the USB port into an adapter that is not higher than 5V 2A, and then plug the adapter into the power port to charge.

### Maintenance

#### Clean the Electrode Surface

- Wipe the residual hydrogel on the electrode surface with a wet paper towel after usage, then wipe the water stains with a dry paper towel.
- Do not wipe with rough objects to avoid scratching the electrode pads.
- cannot be wiped clean, wipe it with a cotton swab and then wipe with a wet tissue. Do not immerse device in water for cleaning.
- Do not immerse device in water for cleaning.